


BOX for SURGICAL INSTRUMENTS

Code **17V64** - Code **17V65** - Code **17V66**

The new surgical boxes have been designed to ergonomically store the surgical instruments. They are made of light but sturdy plastic materials, totally autoclavable at 134°C (274°F) at 2,1 bars for 5 minutes.

The pull-out tray included in all three boxes is screen printed with clear and intuitive symbols to facilitate the identification of instruments and accessories and clearly visualize the recommended surgical sequence.

The handling of the instruments (drills, stops and connections) is very easy, thanks to the vertical position.

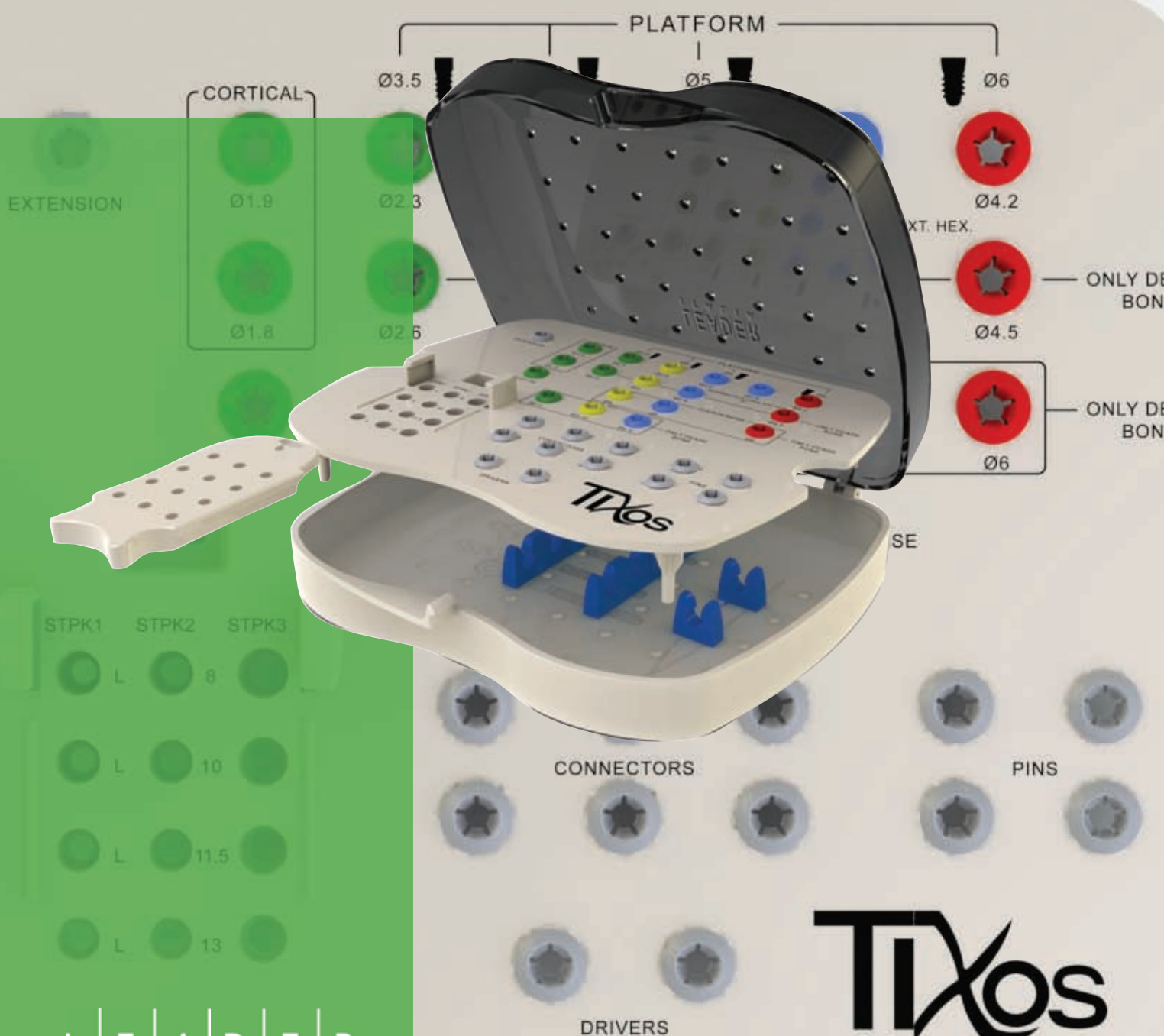

LEADER
ITALIA

TIXOS

LEADER ITALIA srl - Cinisello Balsamo MI - ITALY
ph +39 (0)2 618651 - www.leaderitalia.it - export@leaderitalia.it

BOX for TIXOS implants

Code **17V64** Box for **TIXOS** implants (box without instruments)



Stops housing

The **pull-out tray** included in the boxes shows the colour coding that identifies the available platforms of the implants line.


The screen printed symbols show also the recommended surgical sequence for the four platforms of the implants.

The **great news for the Tixos box** is that the **recommended surgical sequence is different depending on bone quality**. Therefore you will find the indication of the **final drills** to use in **soft bone** and the **extended sequence** with countersinks and bone taps in case of **compact bone**.

A proper housing is specially dedicated to the stops, which can be extracted directly with the drills.

In the **box base**, the complementary surgical accessories can be stored (the ratchet for implants manual insertion, one probe, the tweezers and two open wrenches).

Pull-out tray

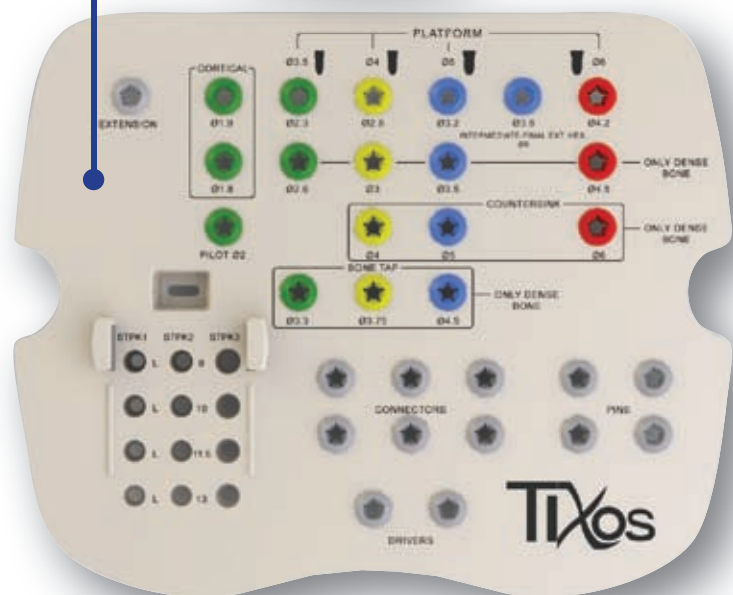

	Platform
Green	Ø 3,5
Yellow	Ø 4,0
Blue	Ø 5,0
Red	Ø 6,0

Box size:
200x165x60 mm

BOX for TIXOS SHORT implants

Code **17V66** Box for **TIXOS SHORT** implants (small box without instruments)



Stops housing

Pull-out tray



As the previous ones, this specific box for Tixos Short is composed by the **pull-out tray** with the screen print of the recommended surgical sequence for the insertion of the Short implant, with a proper housing for the stops which can be extracted directly with the drills.

The **box base** can accommodate the ratchet for the implant manual insertion.

Box size:
110x105x60 mm

BOX for IMPLUS implants

Code **17V65** Box for **IMPLUS** implants (box without instruments)



Stops housing


The **pull-out tray** included in the boxes shows the colour coding that identifies the available platforms of the implants line.

The screen printed symbols show also the recommended surgical sequence for the four platforms of the implants.

A proper housing is specially dedicated to the stops, which can be extracted directly with the drills.

In the **box base**, the complementary surgical accessories can be stored (the ratchet for implants manual insertion, one probe, the tweezers and two open wrenches).

Pull-out tray

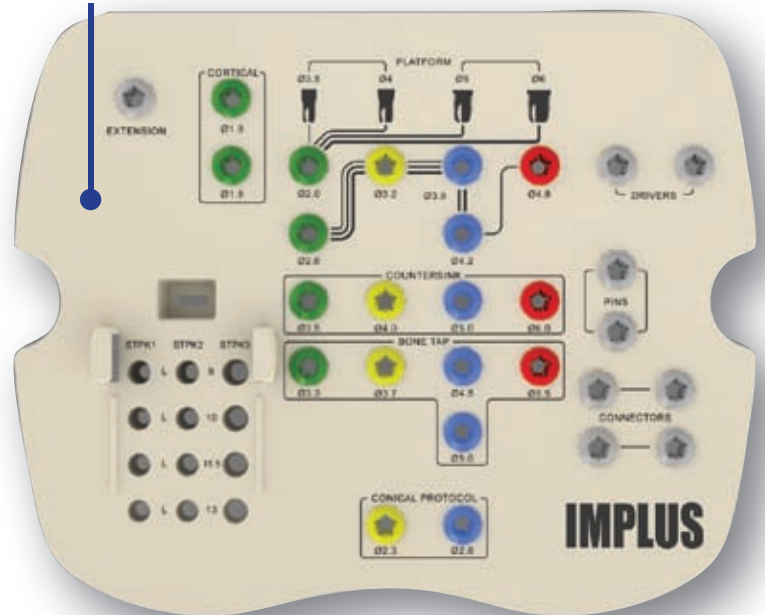

	Platform
Green	ø 3,5
Yellow	ø 4,0
Blue	ø 5,0
Red	ø 6,0

Box size:
200x165x60 mm



Note: the technical suggestions relative to the surgical protocol reported on the boxes can under no circumstances substitute the clinical evaluation of the dentist.